

## Specifications

**Dimensions:** Full size, from 5'5" tall, folded in bag, L x W x H 34" x 12" x 10" Unfolded 38" x 23,5" height variable

Junior; up to 5'5" tall, in bag 3" shorter. Unfolded 3" shorter and 3" narrower

**Shopper** (lowest weight version) - slightly reduced dimensions due to smaller wheels. In the Junior version it is perfect for children (6+) and short-statured people

**Total weight** 35 lbs. w/o battery and seat 24 lbs Shopper 26 lbs (18 lbs)

**Choice of 3 li-Ion batteries,**

Watt hours (Wh)	Weight	Distance on level and 220 lbs load
158 Wh*	2,6 lbs.	9 mls
274,4 Wh*	4 lbs.	15 mls
420 Wh**	5,8 lbs.	22 mls

\* Approved for air travel \*\* Not approved for air travel

**Charger** 100 to 240 V Charging time totally discharged battery 5 / 8/ 12 hours

**Max. recommended load** 320 lbs. plus luggage

**Motor:** Motor: 24 VDC 200 Watt, incorporated in left rear wheel, forward/reverse

**Climbing power** with 220 lbs. load 15 % (Shopper 6%)

**Speed:** 4 mph\* Shopper 2,8 mph


\*Within the EU the top speed for non-licensed motor vehicles operated on public property is limited to 6 kph. Higher speed setting available on request.

**Wheels:** front 8" x 2" rear 7,25" x 2,5" all filled with dense foam

**Brakes:** Each rear wheel with band brake and parking brake

Included in each shipment:

1 TravelScoot with battery of choice and charger, duffel bag, fabric "caddy" pouch (for light items; easy to remove and install), tool kit. Cane holder and basket rail are standard.

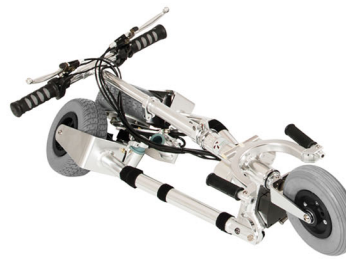

At a sensationally low weight of 35 lbs the **Travelscoot** is worldwide by far the lightest and most compact electric mobility scooter. Ideal for everyone who's feet are not playing along but who is otherwise still in a fairly good shape.

The innovative folding mechanism allows, for the first time, a 3-wheeled scooter to fold down to the dimensions of a 2-wheeler.

Set-up and break-down takes less than a minute. It fits completely into its own carrying bag and can be stored where no mobility scooter has ever fit before!

The **TravelScoot's** snazzy appearance practically eliminates the perception of being a mobility device for the disabled.

Strangers, if they notice you at all, are more likely to give you a "thumbs up" or "cool ride".

**TravelScoot South Africa**  
**Mobility 4 U P.O. Box 416**  
**Yzerfontein Cape Town 7351**  
**South Africa**  
**+27 22 451 2004**  
**+27 84 669 1883**  
**info@mobility4u.co.za**



# TravelScoot

**Re-Mobilize Your Life!**



[www.travelscoot.com](http://www.travelscoot.com)



## Product description:

Our patented, extremely sturdy and lightweight structure is the result of a revolutionary idea combined with high-strength aircraft-grade aluminum. Even though our stated weight limit is 320 lbs, we have many customers who weigh considerably more, and are very willing to accept the modest reduction in climb and range.

Carrying an additional battery "piggyback" style will more than make up for the reduction in range. →

The **TravelScoot**'s incomparable light weight and small size allows it to be taken along everywhere. Whether it's a trip to the corner store, a bus ride into town, or a trip around the world by plane or cruise ship,

*TravelScoot makes it all possible, even if the feet say no. Shopping is no longer painful!*

If walking is difficult, then standing is usually not so great, either. But standing is a daily necessity, whether in line at the bank, at a counter, or just waiting around. If you are on a **TravelScoot**, you are not standing, you **SITTING** in line. You always have your own seat with you on public transport, and don't have to wait for someone to offer you their seat. Just park your **TravelScoot** perpendicular to the direction of travel, and set the hand brakes.

All single-seat electric scooters with more than two wheels are considered handicap vehicles and permitted wherever walking is the only alternative. This includes buildings as well!

Conveniently, the unrivaled reduced vehicle weight results in a marked increase in range.

A mobility scooter will need to be lifted around several times a day, if only to clear a curb or into a vehicle. In some cases, there is no elevator, so you will need to carry your scooter up a flight of stairs.

In such cases, every pound makes a difference!



# We have a comprehensive assortment of accessories available for your **TravelScoot**.

*The convertible trailer can be used like this ..... or like this.*



*For light weight carry on luggage*



*Cruise ship set*



*Foot Rest Support Plate*



*Foldable basket*



*Travel-Set*



*Auxiliary Battery Tray*



*You can choose between 3 different battery sizes / ranges*



*Suspension seat post*